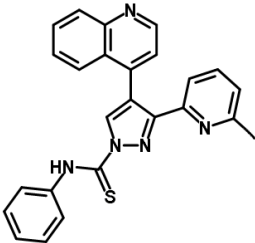




Product Specification Sheet

Product Name	Stemolecule™ A83-01
Description	A83-01 is a selective inhibitor of the transforming growth factor-beta (TGF-β) type I receptor ALK5, the Activin/Nodal receptor ALK4, and the nodal receptor ALK7 ¹ . This molecule is more potent than SB431542 in its inhibition of ALK4, 5, and 7, and only weakly inhibits ALK1, 2, 3, and 6. A83-01 inhibits the TGF-β-induced epithelial-to-mesenchymal transition (EMT) via the inhibition of Smad2 phosphorylation ² .
Catalog Number	04-0014-10
Size	10 mg
Chemical Name	3-(6-Methyl-2-pyridinyl)- <i>N</i> -phenyl-4-(4-quinolinyl)-1 <i>H</i> -pyrazole-1-carbothioamide
Chemical Formula	C ₂₅ H ₁₉ N ₅ S
Structure	
Molecular Weight	421.52
CAS Number	909910-43-6
Purity	Greater than 98% by HPLC analysis
Formulation	Pale yellow solid
Solubility	For a 10 mM concentrated stock solution of A83-01, reconstitute the compound by adding 2.37 ml of DMSO to the entire contents of the vial. If precipitate is observed, warm the solution to 37°C for 2 to 5 minutes. For use in cell culture, warm the medium just prior to adding the reconstituted compound. Once the compound is added, mix and filter-sterilize the medium using a 0.2 μM low-protein binding filter. A83-01 is soluble in DMSO at 50 mM.
Storage and Stability	Store powder at 4°C protected from light. Following reconstitution, store aliquots at -20°C. Stock solutions are stable for 6 months when stored as directed.
Quality Control	The purity of A83-01 was determined by HPLC analysis. The accurate mass was determined by mass spectrometry. No acute cytotoxicity was observed in mouse embryonic stem cells following a 6 hour exposure to 1 nM - 100 μM of A83-01.
References	<ol style="list-style-type: none">Tojo, M., Hamashima, Y., Hanyu, A., Kajimoto, T., Saitoh, M., Miyazono, K., Node, M., and Imamura, T. (2005) The ALK-5 inhibitor A-83-01 inhibits Smad signaling and epithelial-to-mesenchymal transition by transforming growth factor-beta. <i>Cancer Sci</i> 96: 791-800.Li, W., Wei, W., Zhu, S., Zhu, J., Shi, Y., Lin, T., Hao, E., Hayek, A., Deng, H., and Ding, S. (2009) Generation of rat and human induced pluripotent stem cells by combining genetic reprogramming and chemical inhibitors. <i>Cell Stem Cell</i> 4: 16-19.

For research use only. Not for use in diagnostic procedures.

©Stemgent 2012. Unless otherwise noted, Stemgent, Stemgent Logo and all other trademarks are the property of Stemgent, Inc.

Stemgent, Inc | 51 Moulton Street, Cambridge, MA 02138 | www.stemgent.com